





DEFIGARD 5000

A new vision on hospital defibrillation

The DEFIGARD 5000 combines the efficacy of a pulsed biphasic defibrillator (*MULTIPULSE BIOWAVE*[®]) with precision monitoring on the largest colour screen available on the world market.

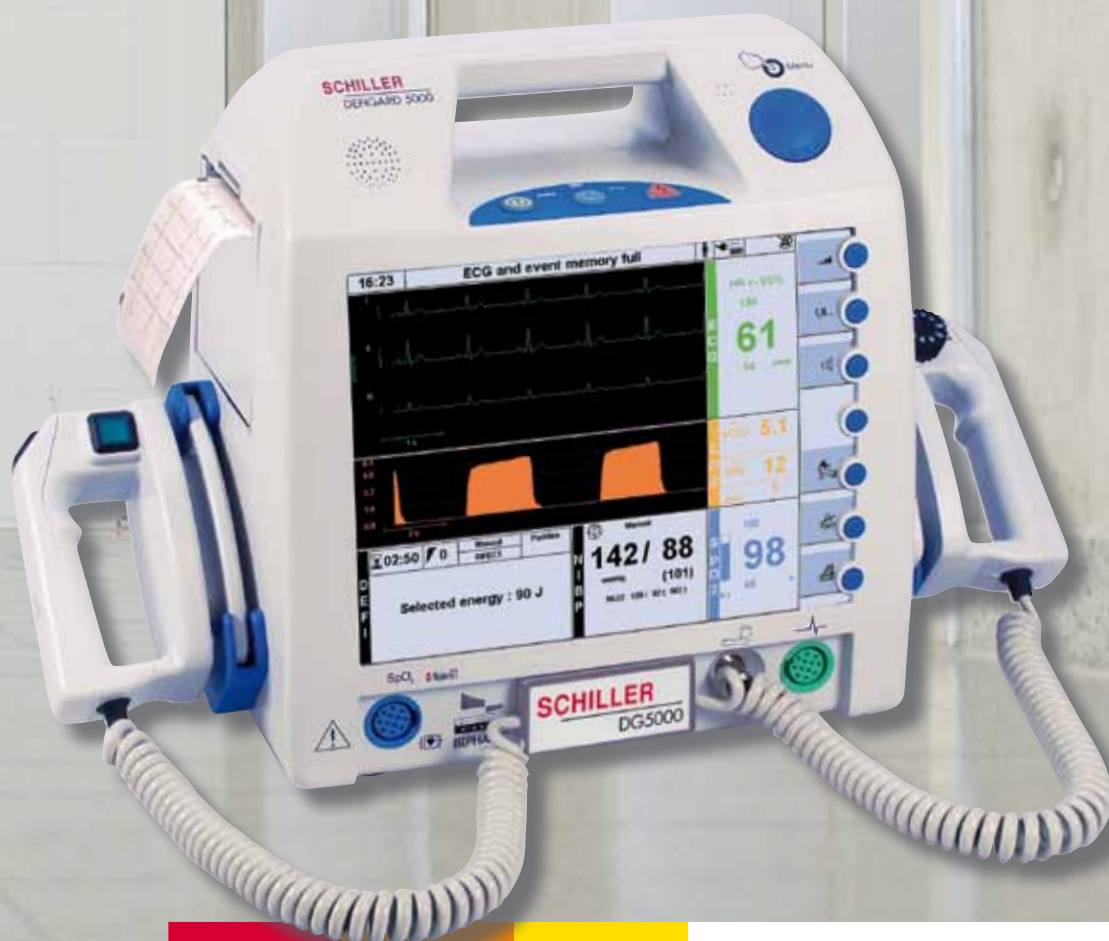

MULTIPULSE
BIOWAVE

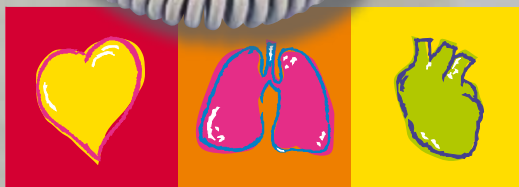

SCHILLER

The Art of Diagnostics

DEFIGARD 5000

Monitor and defibrillator

The large colour screen (10.4") allows the simultaneous display of up to 12 ECG leads. The navigation button and direct soft keys give you quick access to all parameters and menus and make the device extremely user-friendly and easy to use. Thanks to the AED (automated external defibrillator) and manual defibrillation mode, the DEFIGARD 5000 is the ideal device in the case of a cardiac arrest:

- Use of pads, paddles or spoons
- External pacemaker and monitoring functions such as ECG, SpO₂ (Masimo or Nellcor), CO₂ (IRMA) and NIBP make the continuous care and monitoring of the patient possible.



Various sophisticated features available

Defibrillation using the low-energy *MULTIPULSE BIOWAVE*[®] technology. Cassettes permit the fast switching between defibrillation modes

- AED mode with pads
Configurable to 2 to 200 J (adults) or 2 to 70 J (paediatric); automatic detection of paediatric electrodes and energy adjustment
The user is guided through the defibrillation process by voice and text messages
- Manual mode with paddles or pads, energy selection between 2 and 200 J
- Internal defibrillation: 2 to 30 J

DEFIGARD 5000

Pacemaker (option):

- Operational modes Fix, Demand or Overdrive (frequency x3)
- Settings are saved (current and frequency) so that stimulation can be resumed if necessary.

Monitor:

- Large colour LCD (10.4")
- Display: AED instructions, measurement and curves for SpO₂, CO₂, NIBP and up to 12 ECG leads

Printer (option):

- 3-channel printer
- Automatic or manual printout of curves, measurements, events, trends and alarms

Data recording:

- Internal memory allowing storage of 45 minutes of ECG data and 500 events in AED mode, as well as 24 hour graphical and tabular trend data for all parameters
- Download of AED data via USB port

Transmission:

- Transmission of a 12-lead ECG via GSM or a standard modem (optional)
- Ethernet or USB connection for software update



Simultaneous display of 12 ECG leads



Plethysmogram



Asia
SCHILLER Asia-Pacific / Malaysia
52200 Kuala Lumpur, Malaysia
Phone +603 6272 3033
Fax +603 6272 2030
sales@schiller.com.my
www.schiller-asia.com



Austria
SCHILLER Handelsgesellschaft m.b.H.
A-4040 Linz
Phone +43 732 709 90
Fax +43 732 757 000
sales@schiller.at
www.schiller.at



China
SCHILLER Int'l Trading (Shanghai) Co. Ltd.
200335 Shanghai, China
Phone +86-21-62099627
Fax + 86-21-62099623
sales@schiller.com.my
www.schiller.cn



France
SCHILLER Médical S.A.S.
F-67162 Wissembourg/Cedex
Phone +33 3 88 63 36 00
Fax +33 3 88 94 12 82
info@schiller.fr
www.schiller-medical.com



France (distribution France)
SCHILLER France S.A.S.
F-77600 Bussy St Georges
Phone +33 1 64 66 50 00
Fax +33 1 64 66 50 10
infoschiller@schiller-france.fr
www.schiller-france.com



Germany
SCHILLER Medizintechnik GmbH
D-85521 Ottobrunn
Phone +49 89 62 99 81-0
Fax +49 89 609 50 90
info@schillermed.de
www.schillermed.de



Hungary
SCHILLER Diamed Ltd.
H-1141 Budapest
Phone +36 (1) 383-4780 / 460-9491
Fax +36 (1) 383-4778
sales@schiller.at
www.schiller-hungary.hu



India
SCHILLER Healthcare India Pvt. Ltd.
Mumbai - 400 001, India
Phone +91 22 6152 3333/ 2920 9141
Fax +91 22 2920 9142
sales@schillerindia.com
www.schillerindia.com



Japan
SCHILLER Japan, Ltd.
Hiroshima 734-8551
Phone +81 82 250 2055
Fax +81 82 253 1713
koji.maekawa@schiller.jp
www.schiller.jp



Poland
SCHILLER Poland Sp. z o.o.
PL-02-729 Warszawa
Phone +48 22 8432089
Fax +48 22 8432089
schiller@schiller.pl
www.schiller.pl



Russia & C.I.S.
SCHILLER AG Rep. office
125124 Moscow, Russia
Phone +7 (495) 970 11 33
Fax +7 (495) 970 11 33
mail@schiller-ag.com
www.schiller-cis.com



Spain
SCHILLER ESPAÑA, S.A.
E-28230-Las Rozas/Madrid
Phone +34 91 713 01 76
Fax +34 91 355 79 33
schiller@schiller.es
www.schiller.es




Switzerland
SCHILLER-Reomed AG
CH-8953 Dietikon
Phone +41 44 744 30 00
Fax +41 44 740 37 10
sales@schiller-reomed.ch
www.schiller-reomed.ch



Turkey
SCHILLER TÜRKİYE
Okmeydani-Sisli - Istanbul
Phone +90 212 210 8681 (pbx)
Fax +90 212 210 8684
sales@schiller-turkiye.com
www.schiller-turkiye.com



USA
SCHILLER America Inc.
Doral, Florida 33122
Phone +1 786 845 0620
Fax +1 786 845 06 02
sales@schilleramerica.com
www.schilleramerica.com

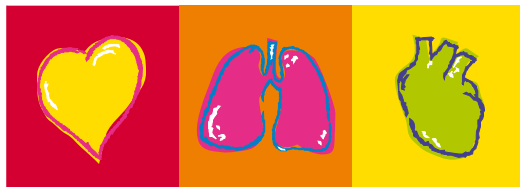

SCHILLER

The Art of Diagnostics

